



Children's
Medical
Ventures

Bendy™ Bumper


The Bendable Positioning Aid
Item 92065

INSTRUCTIONS FOR USE


The Bendy™ Bumper is a flexible, moldable positioning aid that can be shaped and reshaped to any configuration needed by the baby.

1. Remove the outer plastic packaging.

2. Shape the Bendy™ Bumper into a "nest" that will meet the baby's needs. We have suggested a few of the many possible positions.

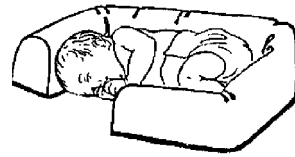

3. Place the baby inside the "nest" in a prone, supine or side-lying position. Consult your hospital physical or occupational therapist for assistance with correct positioning and flexion.

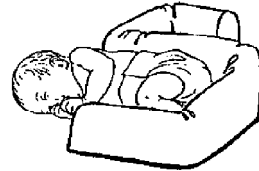


4. If desired, place a blanket over the lower half of the baby and tuck it under the Bendy™ Bumper. This offers the baby more security, prevents the feet from kicking out and helps to maintain proper positioning.

5. To replace the cover: Unfold the open end of the cover. Grasp the foam body, making sure to keep the urethane wrap intact and the inner rod in place. Slide off the cover. Open the end of a clean cover and slide the original foam body with rod into the new cover. Fold ends neatly. **Foam body is single patient use only.** Cloth covers may be laundered and reused according to CDC guidelines.



6. To reuse the rod: Straighten out the Bendy™ Bumper. Pull out the rod (this may be difficult initially due to bends in the rod). Remove the foam body and discard. Place the encapsulated rod on a flat surface and "roll" it back and forth (like a rolling pin) in order to straighten it out as much as possible. Wipe down the tubing with an antibacterial disinfectant. Take a new foam body with cover. Unfold the open end and insert the rod (label side up) into the predrilled hole in the foam. Do not force rod into foam; let it "fall" into the predrilled hole. If it does not go in easily, try rolling the rod again until it is straightened. Hold foam flap up, fold over cover ends neatly.

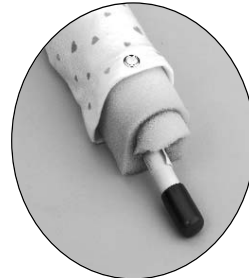

...continued on reverse

HELPING YOU HELP YOUR BABY™



Children's
Medical
Ventures

Bendy™ Bumper


The Bendable Positioning Aid
Item 92065

INSTRUCTIONS FOR USE


The Bendy™ Bumper is a flexible, moldable positioning aid that can be shaped and reshaped to any configuration needed by the baby.

1. Remove the outer plastic packaging.

2. Shape the Bendy™ Bumper into a "nest" that will meet the baby's needs. We have suggested a few of the many possible positions.

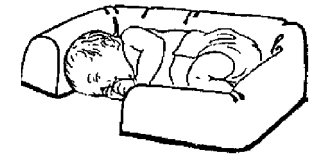

3. Place the baby inside the "nest" in a prone, supine or side-lying position. Consult your hospital physical or occupational therapist for assistance with correct positioning and flexion.

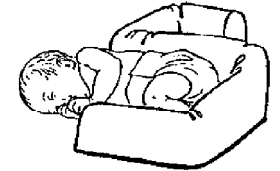


4. If desired, place a blanket over the lower half of the baby and tuck it under the Bendy™ Bumper. This offers the baby more security, prevents the feet from kicking out and helps to maintain proper positioning.

5. To replace the cover: Unfold the open end of the cover. Grasp the foam body, making sure to keep the urethane wrap intact and the inner rod in place. Slide off the cover. Open the end of a clean cover and slide the original foam body with rod into the new cover. Fold ends neatly. **Foam body is single patient use only.** Cloth covers may be laundered and reused according to CDC guidelines.



6. To reuse the rod: Straighten out the Bendy™ Bumper. Pull out the rod (this may be difficult initially due to bends in the rod). Remove the foam body and discard. Place the encapsulated rod on a flat surface and "roll" it back and forth (like a rolling pin) in order to straighten it out as much as possible. Wipe down the tubing with an antibacterial disinfectant. Take a new foam body with cover. Unfold the open end and insert the rod (label side up) into the predrilled hole in the foam. Do not force rod into foam; let it "fall" into the predrilled hole. If it does not go in easily, try rolling the rod again until it is straightened. Hold foam flap up, fold over cover ends neatly.

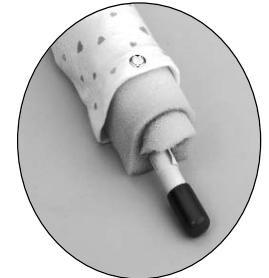

...continued on reverse

HELPING YOU HELP YOUR BABY™



Bendy™ Bumper

INSTRUCTIONS FOR USE...continued

CAUTIONS

- THIS IS NOT A TOY. Use under the direction of a health care professional only.
- Do not place Bendy™ Bumper directly in front of the baby's face or nose or in any location that the baby's airways would be restricted.
- Always use with a cover. If the baby is placed in direct contact with the urethane wrap, it may cause restriction of the airway.
- Do not place on top of baby, lines or tubes.
- The encapsulated rod contains lead. Do not open or cut plastic tubing (see rod label).
- If encapsulated rod cover is damaged in any way, discontinue use and call 1-800-345-6443 for disposal instructions.
- Encapsulated rod may be recycled. DO NOT INCINERATE.

CARE and MAINTENANCE

- **Fabric cover:** Cover may be removed for laundering (see label). After laundering, cover should be stretched lengthwise before replacing.
- **Disposable non-woven cover:** Single patient use only. Discard after use.
- **Foam body:** Single patient use only. Urethane wrap may be wiped down with disinfectant, if necessary (see label). Discard foam body and sleeve between patients or if foam body is heavily soiled. Do not launder or wash foam.
- **Encapsulated rod:** Reusable. Rod must be wiped down with disinfectant between patients (see label).



Respironics Inc.
1001 Murry Ridge Lane, Murrysville, PA 15668
USA 800-345-6443 • Outside USA 1-724-387-4000

www.childmed.com • www.respironics.com

92065-0012
HM 4/24/03

HELPING YOU HELP YOUR BABY™



Bendy™ Bumper

INSTRUCTIONS FOR USE...continued

CAUTIONS

- THIS IS NOT A TOY. Use under the direction of a health care professional only.
- Do not place Bendy™ Bumper directly in front of the baby's face or nose or in any location that the baby's airways would be restricted.
- Always use with a cover. If the baby is placed in direct contact with the urethane wrap, it may cause restriction of the airway.
- Do not place on top of baby, lines or tubes.
- The encapsulated rod contains lead. Do not open or cut plastic tubing (see rod label).
- If encapsulated rod cover is damaged in any way, discontinue use and call 1-800-345-6443 for disposal instructions.
- Encapsulated rod may be recycled. DO NOT INCINERATE.

CARE and MAINTENANCE

- **Fabric cover:** Cover may be removed for laundering (see label). After laundering, cover should be stretched lengthwise before replacing.
- **Disposable non-woven cover:** Single patient use only. Discard after use.
- **Foam body:** Single patient use only. Urethane wrap may be wiped down with disinfectant, if necessary (see label). Discard foam body and sleeve between patients or if foam body is heavily soiled. Do not launder or wash foam.
- **Encapsulated rod:** Reusable. Rod must be wiped down with disinfectant between patients (see label).



Respironics Inc.
1001 Murry Ridge Lane, Murrysville, PA 15668
USA 800-345-6443 • Outside USA 1-724-387-4000

www.childmed.com • www.respironics.com

92065-0012
HM 4/24/03

HELPING YOU HELP YOUR BABY™